

2019 NIA IRP Summer Intern Lab Safety Manual

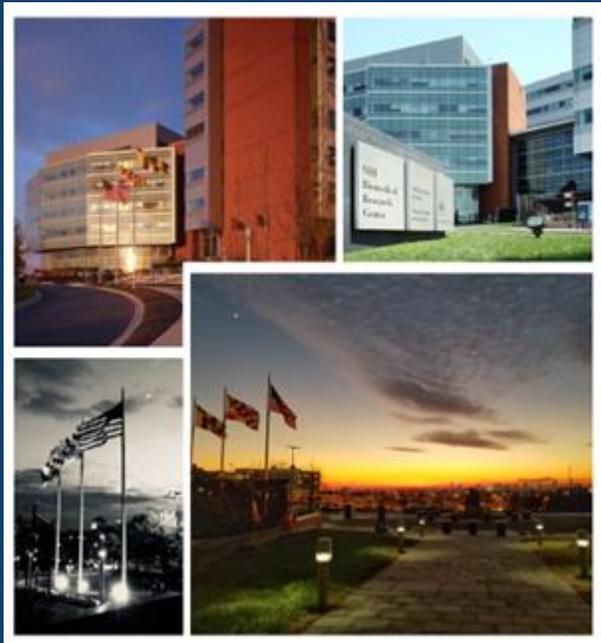


Table of Contents

I. Responsibilities	Page 2
II. Appropriate Dress and Behavior	Page 3
III. Safety Training	Page 4
IV. General Safety Rules	Page 5
V. Working with Chemicals	Page 6
VI. Biological Material	Page 11
VII. Engineering Safety Controls	Page 12
VIII. Personal Protective Equipment	Page 14
IX. Fire Safety	Page 16
X. Ergonomics	Page 17
XI. Preventing Musculoskeletal Disorders	Page 18
XII. Work Related Injuries	Page 19
XIII. NIA IRP Safety Resources	Page 20



NIA IRP ANNUAL SAFETY REFRESHER

Responsibilities

As a summer intern you will work under the supervision of your mentor and various other investigators. However,

YOU ARE ULTIMATELY RESPONSIBLE FOR YOUR SAFETY!

While your duties will vary according to your interest and National Institute on Aging (NIA) IRP Branch assignment, you must follow all safety procedures developed by the National Institutes of Health (NIH) and NIA IRP. All of the policies and procedures are mandated by the federal government through the NIH, the Occupational Safety and Health Administration (OSHA), the Environmental Protection Agency (EPA) and various other agencies.

A working understanding of laboratory safety is integral to any scientific discipline. Aside from injury, loss of resources, and recovery time, laboratory incidents have legal ramifications. An organization can be fined for safety violations even when there are no accidents. Safety knowledge is a skill set that employers expect scientists to have.

Common Causes of Lab Accidents

- Lack of working understanding of hazards
- Improper or unintended use of equipment or reagents
- Inexperienced scientist
- Distractions, lack/loss of attention to task
- Broken, damaged glassware or equipment

You are required to report all work-related accidents, injuries and illnesses to your supervisor, mentor, and the NIA IRP Safety Office. You have the right to know the dangers of the material and equipment used. In addition you are required to wear the Personal Protective Equipment (PPE) provided for your protection. While working in the lab you are expected to dress and behave appropriately. Paying attention to the task you are performing is critical in preventing accidents in the laboratory. Try to do one thing at a time, and do it well. Anticipate and minimize distractions such as phone calls, radios, talking while performing your experiments.

You will be tested on the material in this manual.

Appropriate Dress and Behavior

APPROPRIATE DRESS HELPS TO ENSURE YOUR SAFETY

- Do not wear shorts, sandals, open toed shoes or shoes made of woven material. Sturdy shoes that cover the entire foot are necessary for protection against corrosive materials and mishaps with equipment.
- Protect your skin. Wearing bare midriiffs or low cut pants may result in severe scaring if you or a coworker has an accidental spill or splash.
- Be aware of dangling jewelry, loose or baggy clothing, or long hair that might get caught in equipment.
- If dressed inappropriately, you may be asked to leave the lab and return in proper attire.

Professional standards of personal behavior are expected in any laboratory.

- Avoid distracting or startling other workers.
- Conduct yourself in a responsible manner at all times in the laboratory. Do not engage in practical jokes or horseplay.
- Do not sit on the laboratory benches or tables.
- Do not allow visitors in laboratories where hazardous materials are stored or are in use. Access to the labs should be limited to authorized personnel only.
- Do not eat, drink, chew gum or apply cosmetics in the laboratory.
- Contact lenses are not recommended for laboratory work since that may trap chemicals causing injury to the eye.



Safety Training

Detailed safety and chemical handling procedures will be explained to you in a number of safety training classes, as well as by your mentor/supervisor. You and your mentor will receive a detailed training schedule based on your work assignment.

- **NIA IRP Baltimore On-site Orientation** – *Required for all summer interns:* This information will be presented during the initial orientation. If you begin working before or after the orientation date, a one on one training session may be provided by the NIA IRP Safety Office. Upon completion, you will be given a short quiz on the material. NOTE: As a returning student, you may be instructed to complete NIA IRP Annual On-site Safety Training instead of On-site Orientation.
- **Introduction to Laboratory Safety** – *Required if you are assigned to work in a laboratory.*
- **NIH Laboratory Safety** – *Required if you are assigned to work in a laboratory.*
- **S.T.A.R.S. Learn By doing Lab Safety Training Classroom** – *Required in place of NIH Laboratory Safety if you are 21 or under and assigned to work in a laboratory.*
- **Working Safely with HIV and BBP** – *Required if you are assigned to work in Biosafety Level (BSL) 2 or above lab space, or if you are working with any of the following: human subjects, human blood/body fluids, or human cell lines, non-human primate tissue, select agents, known pathogens and/or any potentially infectious materials.*
- **NIA IRP Annual On-site Safety Training** – *Required every year for returning employees, interns, students, and contractors regardless of job responsibilities.*

All new summer interns will be required to complete the On-site Orientation. Depending upon your job responsibilities you may also be required to take Introduction to Laboratory Safety, NIH Laboratory Safety, S.T.A.R.S., and Working Safely with HIV and Other Bloodborne Pathogens. Returning summer interns may be required to complete the annual safety refresher titled “NIA IRP Annual On-site Safety Training.” Additional training may be required if your responsibilities have changed.

General Safety Rules

- Become familiar with the locations and use of all safety equipment including eye wash stations, safety showers, spill control materials, fire blankets, fire extinguishers, fire alarms and fire exits.
- Keep the aisles clear.
- Do not work alone.
- Observe good housekeeping practices. Clean up spills immediately.
- Examine glassware before each use. Do not use dirty, chipped, broken or cracked glassware in any procedure.
- Do not prepare or store food or drink in any laboratory.
- Do not taste, touch or smell any reagents.
- Wear appropriate gloves and eye protection while handling infectious or toxic materials and animals.
- Wear appropriate clothing to cover exposure skin. Open-toed shoes, sandals, shorts and other apparel that leave skin exposed are not appropriate for wear in any laboratory.
- Wear laboratory coats and other personal protective equipment/clothing (PPE) while performing laboratory activities. **Do not wear PPE** outside of the laboratory.
- Never pipet by mouth.
- Always wash your hands before leaving the laboratory. This will prevent contamination of your food or person.
- Remove gloves carefully (discard disposable gloves immediately after use); thoroughly wash hands and forearms upon completion of work and before leaving the laboratory.
- Use 10% bleach solution (or equivalent) prepared daily to disinfect work station at beginning and end of work.

Working with Chemicals

There are hazards associated with most chemicals. During your safety trainings at the NIA IRP, you will be given information on hazard identification, safety protocols and emergency procedures. Please remember even though you will be given training, personal protective equipment and supervision, *you are responsible for your safety*. In addition to the general safety protocols listed above, there are additional safety measures needed when working with chemicals:

- Do not store chemicals near heat sources such as ovens or steam pipes. Also, do not store chemicals in direct sunlight.
- Do not use lab benches as permanent storage for chemicals. Chemicals should be stored where they cannot be easily knocked over, incompatible chemicals are not stored alongside one another, and the chemicals are protected in the event of a fire.
- Each chemical must have a proper designated storage location and be returned there after use.
- Inspect your chemical containers routinely for any signs of deterioration and for the integrity of the label. State and federal law requires that **all** chemicals must be clearly labeled.
- Do not store any chemicals containers directly on the floor; where appropriate, use secondary containment.
- Do not use a Chemical Fume Hood (CFH) as a permanent storage location for chemicals, with the exception of particularly odorous chemicals that may require ventilation. The more containers, boxes, equipment, and other items that are stored in a CFH, the greater the likelihood of having chemical vapors being drawn back into the room.
- Chemicals that require refrigeration must be sealed with tight-fitting caps and kept in lab safe refrigerators. Lab safe refrigerators/freezers must be used for flammables requiring cold storage.
- Do not store chemicals above eye level. If the container breaks, the contents can fall onto your face and upper body.
- Use a chemical fume hood when opening, pouring, or handling hazardous chemicals. NOTE: do not use biological safety cabinet (BSC) for chemical use.
- Transport laboratory chemicals using bottle carriers and suitable carts.
- The identification and disposal of unlabeled chemical containers is very expensive. Label all chemical containers/solutions in the laboratory with the following information:
 - The name of the chemical or stock solution
 - The date of preparation
 - Concentration
 - Your initials
- Ensure unimpeded access to safety showers and eyewash stations. Test and flush eyewash stations weekly.

- Do not dispose of chemicals down the drain or by evaporation.
- Properly collect, tag, and date waste. Keep chemical waste containers closed/sealed. Use drip pans under waste collection containers to prevent spills.
- Additional information on specific chemical handling, storage and use may be found on chemical labels, Safety Data Sheets (available online and lab binders) and the safety office.

Flammable/Combustible Liquids

Flammable liquids are indeed the most common chemicals found in a laboratory. The primary hazard associated with flammable liquids is, of course, their ability to readily ignite and burn. One should note that it is the vapor of a flammable liquid, not the liquid itself that ignites and causes a fire. The rate at which a liquid vaporizes is a function of its vapor pressure. Vapor pressure increases rapidly as the temperature is raised as does the evaporation rate. A reduced-pressure environment also accelerates the rate of evaporation. The flash point of a liquid is the lowest temperature at which a liquid gives off enough vapor to form a mixture that will burst into flame if given a source of ignition.

- Flammable liquids that are not in active use must be stored in safe containers inside fire resistant storage cabinets designed for flammables.
- Use flammables only in areas free of ignition sources.
- Do not heat flammables with an open flame. Instead, use steam baths, water baths, oil baths, hot air baths, sand baths or heating mantles.
- Never store flammable chemicals in a standard household refrigerator. There are several ignition sources located inside a standard refrigerator that can set off a fire or violent explosion. Flammables can only be stored cold in a lab safe or explosion-proof refrigerator. Another alternative is to use an ice bath to chill the chemicals.

Oxidizers

Oxidizers or oxidizing agents present fire and explosion hazards on contact with combustible materials. Depending on the class, an oxidizing material may increase

the burning rate of combustibles with which it comes in contact; cause the spontaneous ignition of combustibles with which it comes in contact; or undergo an explosive reaction when exposed to heat, shock, or friction. Oxidizers are generally corrosive.

- Store oxidizers away from flammables, organic compounds, and combustible materials.
- Strong oxidizing agents like chromic acid should be stored in glass or some other inert container, preferably unbreakable. Corks and rubber stoppers should not be used.

- Reaction vessels containing appreciable amounts of oxidizing material should never be heated in oil baths, but rather on a heating mantle or sand bath.

Corrosives

Corrosives are most commonly acids and alkalis, but many other materials can be severely damaging to living tissue. Corrosives can cause visible destruction or irreversible alterations at the site of contact. Inhalation of the vapor or mist can cause severe bronchial irritation. Corrosives are particularly damaging to the skin and eyes.

Please be aware that certain substances considered non-corrosive in their natural dry state are corrosive when wet such as when in contact with moist skin or mucus membranes. Examples of these materials are lithium chloride, halogen fluorides, and allyl iodide. Sulfuric acid is a very strong dehydrating agent and nitric acid is a strong oxidizing agent. Dehydrating agents can cause severe burns to the eyes due to their affinity for water.

- Always store acids separately from bases. Also, store acids in acid storage cabinets away from flammables since many acids are also strong oxidizers.
- Do not work with corrosives unless an emergency shower and continuous flow eyewash are available. Know the location and correct operation of these essential pieces of safety equipment.
- Add acid to water, but never add water to acid. This is to prevent splashing from the acid due to the generation of excessive heat as the two substances mix.
- Never store corrosives above eye level. Store on a low shelf or cabinet.
- It is a good practice to store corrosives in a tray or bucket to contain any leakage.

Water Reactive Materials

When water reactive materials come in contact with water, one or more of the following can occur: liberation of heat which may cause ignition of the chemical itself if it is flammable, or ignition of flammables that are stored nearby; release of a flammable, toxic, or strong oxidizing gas; release of metal oxide fumes; and formation of corrosive acids. Water reactive chemicals can be particularly hazardous to firefighting personnel responding to a fire in a lab, because water is the most commonly used fire extinguishing medium.

Peroxide-Forming Materials

Peroxides are very unstable and some chemicals that can form them are commonly used in laboratories. This makes peroxide-forming materials some of the most hazardous substances found in a lab. Peroxide-forming materials are chemicals that react with air, moisture, or impurities to form peroxides. The tendency to form peroxides by most of these materials is greatly increased by evaporation or distillation. Organic peroxides are extremely sensitive to shock, sparks, heat, friction, impact, and light. Many peroxides formed from materials used in laboratories are more shock sensitive than TNT. Just the friction from unscrewing the cap of a container of ether that

has peroxides in it can provide enough energy to cause a severe explosion. Examples of peroxide-forming materials: diisopropyl ether, divinylacetylene sodium amide, potassium amide, dioxane, diethyl ether, tetrahydrofuran vinyl ethers, butadiene, vinylpyridine, acrylonitrile, styrene.

- Do not open the chemical container if peroxide formation is suspected. The act of opening the container could be sufficient to cause a severe explosion. If you notice crystal formation in the container or around the cap, do not attempt to open or move the container. Call the Safety Office or inform your mentor.
- Visually inspect liquid peroxide-forming materials for crystals or unusual viscosity before opening. Pay special attention to the area around the cap. Peroxides usually form upon evaporation, so they will most likely be formed on the threads under the cap.
- Date all peroxide forming materials with the date received, and the expected shelf life. Chemicals such as diisopropyl ether, divinyl acetylene, sodium amide, and vinylidene chloride should be discarded after three months.
- Store all peroxide-forming materials away from heat, sunlight, and sources of ignition. Sunlight accelerates the formation of peroxides.
- Secure the lids and caps on these containers to discourage the evaporation and concentration of these chemicals.
- Never store peroxide-forming materials in glass containers with screw cap lids or glass stoppers. Friction and grinding must be avoided. Also, never store these chemicals in a clear glass bottle where they would be exposed to light.

Compressed Gas Cylinders

Cylinders of compressed gases can pose a chemical as well as a physical hazard. If the valve were to break off a cylinder, the amount of force present could propel the cylinder through a brick wall. For example, a cylinder of compressed breathing air used by SCUBA divers has the explosive force of 1 ½ pounds of TNT.

- Always use the appropriate regulator on a cylinder. If a regulator will not fit a cylinder's valve, replace the cylinder, not the regulator. Do not attempt to adapt or modify a regulator to fit a cylinder it was not designed for. Regulators are designed to fit only specific cylinder valves to avoid improper use.
- Inspect regulators, pressure relief devices, valves, cylinder connections, and hose lines frequently for damage.
- Never use a cylinder that cannot be positively identified. Color coding is not a reliable way of identifying a cylinder because the colors can vary from supplier to supplier.
- Do not use oil or grease on any cylinder component of an oxidizing gas because a fire or explosion can result.
- Never transfer gases from one cylinder to another. The gas may be incompatible with the residual gas remaining in the cylinder or may be incompatible with the cylinder material.
- Place all cylinders so that the main valve is always accessible.

- Close the main cylinder valve whenever the cylinder is not in use.
- Remove regulators from unused cylinders and always put the safety cap in place to protect the valve.
- Always secure cylinders, whether empty or full, to prevent them from falling over and damaging the valve (or falling on your foot). Secure cylinders by chaining or strapping them to a wall, lab bench, or other fixed support.
- To transport a cylinder, put on the safety cap and strap the cylinder to a hand truck in an upright position. Never roll a cylinder.
- Always clearly mark empty cylinders and store them separately.
- Be careful while handling compressed gas cylinders and never drop or strike a cylinder against anything.
- Use only wrenches or other tools supplied by the cylinder supplier to open a valve.
- Open cylinder valves slowly.

Globally Harmonized System		
 <p>Expanding Bomb This chemical can blow up.</p>	 <p>Flame Flammable chemicals can catch fire easily and burst into flames.</p>	 <p>Environment This chemical can kill fish and other life that live in water.</p>
 <p>Corrosive Corrosive chemicals can cause serious damage to the skin and eyes. They can also eat away clothing, metal and work surfaces.</p>	 <p>Gas Cylinder This chemical can explode, rocket or harm health if the cylinder is heated, ruptured or leaking.</p>	 <p>Flame Over Circle Oxidizing chemicals can react with other materials causing them to burn or explode.</p>
 <p>Health Hazard Prolonged exposure this chemical may cause health problems such as cancer or birth defects. Some chemicals may cause asthma or damage to specific organs of the body.</p>	 <p>Exclamation Mark This chemical may cause immediate health effects such as skin rashes or respiratory irritation. Some chemicals may damage the ozone layer.</p>	 <p>Skull & Crossbones Exposure to this chemical can cause immediate and possibly serious health problems.</p>

Biological Hazards

Working with bio-hazardous material requires special safety precautions. Viruses, bacteria and other microorganisms can enter the body through contact with non-intact skin or mucus membranes, by inhalation, ingestion and inoculation. Protection is provided by primary barriers such as biological safety cabinets, personal protective equipment and special procedures.

Standard Microbiological Practices

- Persons wash their hands after they handle viable materials and animals, after removing gloves, and before leaving the laboratory.
- Eating, drinking, handling contact lenses, and applying cosmetics are not permitted in work areas where there is reasonable likelihood of exposure to potentially infectious materials.
- Persons who wear contact lenses in laboratories should also wear goggles or a face shield.
- Food is stored outside the work area in cabinets or refrigerators labeled and used for this purpose only.
- Mouth pipetting is prohibited; mechanical pipetting devices are used.
- All procedures are performed carefully to minimize the creation of splashes or aerosols.
- Work surfaces are decontaminated at least once a day and after any spill of viable material.
- All cultures, stocks, and other regulated wastes are decontaminated before disposal by an approved decontamination method.
- It is recommended that laboratory coats, gowns, or uniforms be worn to prevent contamination or soiling of street clothes.
- Gloves should be worn if the skin on the hands is broken or if a rash exists.
- Protective eyewear should be worn for anticipated splashes of microorganisms or other hazardous materials to the face.
- Protective laboratory coats, gowns, smocks, or uniforms designated for lab use are worn while in the laboratory. This protective clothing is removed and left in the laboratory before leaving for non-laboratory areas (e.g., cafeteria, library, administrative offices). **DO NOT WEAR PPE OUTSIDE THE LAB.**
- All protective clothing is either disposed of in the laboratory or laundered by the institution; it should never be taken home by personnel.
- Gloves are worn when handling infected animals and when hands may contact infectious materials, contaminated surfaces, or equipment. Wearing two pairs of gloves may be appropriate; if a spill or splatter occurs; the hand will be protected after the contaminated glove is removed.
- Gloves are disposed of when contaminated, removed when work with infectious materials is completed, and are NOT worn outside the laboratory.
- Disposable gloves are not washed or reused.

Engineering Controls

Engineering controls are the first line of defense against chemical and biological hazards. When working with certain hazardous chemicals a chemical fume hood (CFH) is used. When not in use, flammable and combustible materials are kept in a flammable safety cabinet. Protection from potentially infectious material is provided by a biological safety cabinet (BSC).

Chemical Fume Hood (CFH)

Work that involves hazards and noxious materials which are toxic, odoriferous, volatile or harmful shall be conducted within a laboratory hood. The primary purpose of a laboratory hood is to keep toxic or irritating vapors and fumes out of the general laboratory working area. A secondary purpose is to serve as a shield between the worker and equipment being used when there is the possibility of an explosive reaction. This is done by lowering the sash of the hood.

- Always assure the hood is operational before initiating an experiment.
- Do not use infectious material in a chemical fume hood.
- Only items necessary to perform the present experiment should be in the hood. The more equipment in the hood, the greater the air turbulence and the chance for gaseous escape into the lab.
- Substitute toxic chemicals with less hazardous materials whenever possible.
- Keep fume hood exhaust fans on at all times.
- Perform all work six inches inside the hood.
- Never place your head inside the hood.
- Keep the hood sash closed as much as possible at all times to ensure the optimum face velocity and to minimize energy usage.
- Keep lab doors closed to ensure negative room pressure to the corridor and proper air flow into the hood.
- Do not store chemicals in the hood.
- Keep the slots of the baffle free of obstruction.
- Do not use the hood as a waste disposal mechanism (e.g., for evaporation of chemicals).
- Avoid rapid movements in front of the hood including opening and closing the fume hood sash rapidly and swift arm and body movements in front of or inside the hood. These actions may increase turbulence and reduce the effectiveness of fume hood containment.
- Do not override or disable mechanical stops on the sash.

Biological Safety Cabinet (BSC)

A biological safety cabinet (BSC) works by pulling room air across the cabinet work surface and through a high efficiency particulate air (HEPA) filter. Part of the air is re-circulated within the BSC and part is returned to the room. Because of this room air return, do not use hazardous chemicals in a BSC.

- Minimize movement in and out of the BSC. Do not use a sweeping motion. Move arms in and out slowly, perpendicular to the face opening.
- As much as possible, reduce other activities in the room. Minimize rapid movement near the BSC as well as open/closing of the room door
- Materials or equipment placed inside the cabinet may disrupt airflow, resulting in turbulence, possible cross contamination, and/or breach of containment. Extra supplies should be stored outside the cabinet. Only materials and equipment required for the immediate work should be placed in the BSC.
- Avoid blocking the front or back grille. Do not place papers, pipetting devices, etc. on the grille. Do not rest arms on the front grille while working. This may draw room air directly into the cabinet. To prevent this, raise arms slightly above the grille surface.
- Allow the cabinet to stabilize and air sweep to occur. To remove particles in the cabinet, allow the blowers to run at least 15 minutes prior to beginning work. Delay manipulation of material for one minute after placing hands/arms in the BSC to air sweep the arms and hands of surface microbial contaminants.
- Work 4 - 6 inches inside the cabinet.
- An open flame creates air turbulence. Do not use a Bunsen burner. If absolutely necessary, a touch-plate microburner may be used (i.e., Touch-O-Matic).
- To reduce contamination, wipe containers and material with 70% EtOH prior to placing them in the cabinet.
- Work should flow from the clean to the “dirty” areas of the cabinet. Clean cultures on one side can be inoculated in the center and contaminated pipets placed in a container on the other side of the BSC.
- Decontamination of the cabinet interior should occur prior to and upon completion of work.
- Wipe down the interior of the glass, the work surface, cabinet sides and back. Take care not to wet the supply filter diffuser. A HEPA filter loses integrity when wet.

Personal Protective Equipment

Personal Protective equipment (PPE) protects only if it is worn and used properly. In prior years, a chemical splash has been the most common type of student injury. In those cases proper PPE was not worn.

Protective Eyewear

- Goggles provide the best all-around protection against chemical splashes, vapors, dusts, and mists.
- Goggles that have indirect vents or are non-vented provide the most protection.
- Standard safety glasses provide protection against impact.
- If using a laser or UV lights, wear safety glasses or goggles which provide protection against the specific wavelength.
- Remember, prescription glasses do not provide adequate protection in a laboratory setting. Prescription safety glasses can be purchased from most opticians.
- Contact lenses should not be worn in a laboratory because they can trap contaminants behind them and reduce or eliminate the effectiveness of flushing with water from an eyewash station. Contact lenses may also increase the amount of chemicals trapped on the surface of the eye and decrease removal of the chemical by tearing. If it is necessary to wear contact lenses in a lab, wear protective goggles at all times.

Protective Gloves

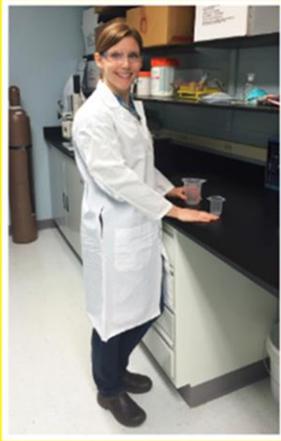
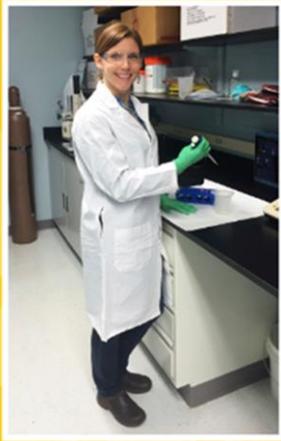
- Any glove can be permeated by chemicals; refer to SDS for best options. The rate at which this occurs depends on the composition of the glove, the chemicals present and their concentration, and the exposure time to the glove. If you are not certain which type of glove provides you with the protection you need, check with your mentor or contact the Safety Office at 410-558-8311.
- If chemical contact occurs or suspected, wash hands and replace gloves immediately.
- Remove gloves and wash hands regularly before using the telephone or opening doors to prevent the spread of contamination.
- Check gloves for cracks, tears, and holes. If compromised, dispose immediately.
- Butyl, neoprene, and nitrile gloves are resistant to most chemicals, e.g., alcohols, aldehydes, ketones, most inorganic acids, and most caustics.
- Disposable latex and vinyl gloves protect against some chemicals, most aqueous solutions, and microorganisms and reduce risk of product contamination.
Disposable gloves are NOT to be reused and must be discarded immediately after use.
- Leather and some knit gloves will protect against cuts, abrasions, and bites/scratches, but not against chemicals.

- Temperature-resistant gloves protect against cryogenic liquids, flames, and high temperatures.

Other Protective Equipment

- The primary purpose of a lab coat is to protect against splashes and spills. A lab coat should be nonflammable, where necessary, and should be easily removed.
- Face shields can protect against impact, dust, particulates, and chemical/biological splashes for the face, eyes, and throat. However, chemical vapors can still travel under and around a face shield. Wear protective eyewear such as goggles underneath a face shield if there is a danger of aerosol formation or splashes that may come under the shield.

Personal Protection Equipment in NIH Laboratories

minimum PPE for ALL LAB WORK	minimum PPE for BSL2 LAB WORK	minimum PPE for BSL2/3 LAB WORK
		
		
<ul style="list-style-type: none"> • Lab coat, plus; • Long pants, plus; • Closed-toed shoes, plus; • Proper eye protection for the task, plus; • Proper gloves for the task. 	<ul style="list-style-type: none"> • Lab coat, long pants & closed-toed shoes, plus; • Safety glasses (minimum) or eye protection appropriate for the task, plus; • Gloves (minimum) or hand protection appropriate for the task. 	<ul style="list-style-type: none"> • Closed-front gown, plus; • Long pants, & closed-toed shoes, plus; • Safety glasses (minimum) or eye protection appropriated for the task, plus; • Double gloves (minimum) or hand protection appropriate for the task.
<p>! = additional PPE appropriate for the task as determined by risk assessment</p>		
<p>NOTE: Do not forget about hazards at adjacent lab benches!</p>		
<p>The National Institutes of Health Office of Research Services Division of Occupational Health and Safety</p>	<p>Reference: NIH Personal Protective Equipment Policy, 2015</p>	

Fire Safety

Most laboratories contain hazardous chemicals which may be flammable. It is important to follow the fire safety guidelines. These will be reviewed in greater detail during your safety training.

- Know the location of the fire alarm pull stations.
- Be familiar with two safe exits from your work place.
- In case of a fire, **DO NOT USE THE ELEVATORS**. Use the stairs to exit.
- Do not use a fire extinguisher unless you have completed Fire Extinguisher Training. Always pull fire alarm before trying to extinguish.
- If your clothes catch fire, do not run. Running will fan the flames. Drop to the ground and begin to roll to smother the flames. While rolling cover your nose and mouth with your hands. This will prevent inhalation of smoke and fire

WHEN THE ALARM SOUNDS...



■ Red Box Fire Alarm

- Everyone evacuates
- Alarm notifies local fire department
- High Rise Evacuation
 - Alarm area, above and below = Evacuate
 - All other areas = Await further instructions

■ Yellow Lab Alarm

- Lab occupants evacuate
- Alarm notifies GSH

**NOTE: Do not use elevators during emergencies
All elevators recall to exit level**

Ergonomics

Ergonomics is the science of fitting jobs to people. It focuses on designing workstations, tools and work tasks for safety, efficiency and comfort. Effective ergonomic design reduces discomfort and injuries and increases job satisfaction and productivity.

Common job injuries that can be reduced with good ergonomics

Injuries affecting muscles, tendons, ligaments, joints, nerves and discs can be reduced or eliminated with ergonomics. These injuries are frequently called Musculoskeletal Disorders (MSD's), Repetitive Strain Injuries (RSI's) and Cumulative Trauma Disorders (CTD's).

Common ergonomic risk factors

Some risk factors have nothing to do with the work environment and can include medical conditions such as rheumatoid arthritis, diabetes, and hormonal imbalances, as well as hobbies involving repetitive motions or awkward postures. Work related risk factors may include:

- Awkward body postures: maintaining an awkward, unsupported or fixed posture such as bending, reaching or twisting for prolonged periods.
- Excessive repetition: doing the same motions over and over again.
- Excessive force: physical exertion or pressure on any part of the body while lifting, pushing, pulling or gripping a tool.
- Contact stress: pressure by tools, edges or hard surfaces on soft tissues of the body. The palms and elbows are often at high risk.
- Vibration: using vibrating or impact tools and equipment.

Common symptoms of Musculoskeletal Disorders

Symptoms that go away overnight are usually a sign of fatigue. Symptoms that don't go away may indicate a more serious problem. Employees who have such symptoms should seek medical advice early. MSD's are easier to treat in the early stages. Ignoring them can result in serious injury. Symptoms frequently include:

- Pain
- Numbness and tingling
- Stiffness or cramping
- Inability to hold objects or loss of strength

Preventing Musculoskeletal Disorders

Prevention:

- Warm up and stretch before starting activities that are repetitive, static or prolonged.
- Take frequent breaks from any sustained posture every 20-30 minutes and stretch stiff muscles.
- Respect pain. Change positions or stop whenever activities cause pain.
- Recognize early signs of the inflammatory process, and treat early.
- Only use splints and wrist supports after instruction by your physician or therapist.

Posture:

- Maintain erect position of back and neck with shoulders relaxed. Minimize twisting and bending motions. Position equipment and work tasks so that your body is directly in front of and close to your major work tasks.
- Use proper positioning during all activities. Keep upper arms close to the body, elbows at 100 degrees, forearms neutral (thumb toward ceiling), and wrist straight. Keep feet flat on the floor when seated by proper adjustment of your chair, or use of a footrest.
- Keep wrists as neutral as possible. Avoid extreme motions. There is a safe zone of movement for your wrist. This zone is about 15 degrees in all directions.
- Avoid bending neck forward for prolonged periods of time. If typing from a manuscript, place the document on a holder beside or below your computer screen.
- Avoid static positions for prolonged periods. Muscles fatigue faster when they are held in one position. Keep moving to increase your blood circulation.

Task Modification:

- Whenever possible, alternate activities frequently throughout the day. Rotate heavy and/or repetitive tasks with lighter, less repetitive tasks.
- If symptoms become worse, or a specific activity consistently causes discomfort, reassess the task setup and look for alternative methods.
- Avoid repetitive or prolonged grip/pre-tension activities. Avoid pinching with wrist flexion or wrist deviation (bending to side). Take frequent breaks to stretch and rest hands.
- Avoid tugging, jerking, or pounding with the hand.

Work Related Injuries

For all work related injuries and exposures, complete the following steps:

1

Immediately initiate first aid at the work site

- Skin – Cleanse and bandage
- Skin (Chemical/Biological Exposure) – Wash with soap and water for 15 minutes
- Eyes and mucous membranes – Irrigate for 15 minutes, water or saline

2

**Notify your supervisor, if he or she is readily available
Notify Safety Office at 410-558-8311**

**Report to Occupational Medical Services (OMS) Baltimore
Monday thru Friday – 8am to 4pm – BRC 1B210
Telephone: 443-740-2309**

- This should be accomplished within one hour of exposure. For emergency transport dial 9-911 (from BRC extension)
- If OMS is closed, during normal business hours promptly called OMS Bethesda at **301-496-4411**; after hours, holidays & weekends call the Clinical Center Operator at **301-496-1211** and ask them to page an OMS physician immediately

3

Division of Occupational Health and Safety



DOHS Baltimore
February 2015



NIA IRP SAFETY RESOURCES



Your lab/area
written policies
and procedures



Specific protocol
safety
precautions



Your safety
committee
representatives



**NIA Safety
Offices – BRC
04C216 &
02B209**



[NIA IRP Safety
Kiosk](#)

9-911 Alerts
24/7 xSAFE (7233)
NIAIRPBaltimoreSafety
Office@mail.nih.gov

NIA IRP SAFETY RESOURCES

DIVISION OF OCCUPATIONAL HEALTH & SAFETY
TEAM BALTIMORE

National Institute on
Aging

Rhonda Walther
Safety Manager

Sarah Ball
Safety Specialist

Brian Farrell
Chemical Hygienist

NIAIRPBaltimoreSafetyOffice@mail.nih.gov
Hotline: 410-558-8300 xSAFE (x7233)



National Institute on
Drug Abuse

Delores Dobson
Lead Safety Manager

Carrie Wertheim
Safety Specialist

Charlie Kowalski
Chemical Hygienist

NIDASafety@mail.nih.gov
Hotline: 443-740-2400



**NIA SAFETY
COMMITTEE**



**NIA SAFETY
OFFICE**